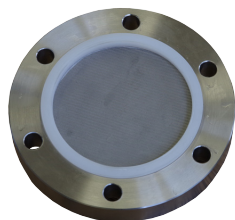


**separchrom**<sup>PC06</sup>

## HPLC COLUMNS FOR PREP CHROMATOGRAPHY

**Columns for Medium Pressure without bed compression**

**Separchrom PC06** are intended for applications where sorbent does not change the volume during the use. Columns are equipped with frits and distributors incorporated in to input and output flanges. They are made of polished stainless steel tubes and equipped with non woven polypropylene frits with porosity 10 µm. Flow distributors are incorporated in to column flanges.

## Columns for laboratory applications

Type of column	I.D. mm	Length mm	Pressure bar	Note
SEPARCHROM PC06 -25/50	25	50	100	
SEPARCHROM PC06 -25/250	25	250	100	
SEPARCHROM PC06 -50/250	50	250	50	
SEPARCHROM PC06 -50/500	50	500	50	
SEPARCHROM PC06 -100/250	100	250	40	
SEPARCHROM PC06 -100/500	100	500	40	

## Columns for industrial applications

Type of column	I.D. mm	Length mm	Pressure bar	Note
SEPARCHROM PC06-100/250	100	250	40	
SEPARCHROM PC06-100/500	100	500	40	
SEPARCHROM PC06-150/300	150	300	40	
SEPARCHROM PC06-150/600	150	600	40	
SEPARCHROM PC06-200/400	200	400	35	
SEPARCHROM PC06-200/800	200	800	35	
SEPARCHROM PC06-300/300	300	300	30	
SEPARCHROM PC06-300/600	300	600	30	
SEPARCHROM PC06-400/200	400	200	20	
SEPARCHROM PC06-400/400	400	400	20	
SEPARCHROM PC06-400/800	400	800	20	
SEPARCHROM PC06-500/500	500	500	10	
SEPARCHROM PC06-500/1000	500	1000	10	
SEPARCHROM PC06-600/600	600	600	5	
SEPARCHROM PC06-600/1200	600	1200	5	